

## ELASTIC LUMBOSACRAL, SACROILIAC SUPPORTS



### MODEL 881-HS

Comfortable construction with a wide range of adjustability. 10" elastic back panel

with 2 pairs of removable steels for rigid support. Velfoam sides and front, tapers from 10" to 7", with front contact closure. Tapered elastic side pulls with hook and loop attachments.

Sizes - Regular: 32"-40" hips  
Large: 42"-48" hips



### MODEL 867

Simple yet effective design provides for comfortable wear as well as economy. 10" back panel with casings for four removable 9" stays if desired (not included). Two-ply elastic sides taper to 6" in front. Two wrap around elastic support straps. Hook and loop fastening and attachment.

Sizes - Small: 30"-34" hips  
Medium: 35"-38" hips  
Large: 39"-42" hips  
X-Large: 43"-46" hips



### MODEL 574

For all figure types; indicated for mild to moderate sacroiliac strain and sprain. White power net elastic, 6" in back, 4" in front, with hook and loop fastening. Removable sacroiliac support pad. Wide elastic encircling pull with hook and loop fastening. Flexible boning in back.

Sizes - Small: 26"-30" hips  
Medium: 32"-36" hips  
Large: 38"-42" hips  
X-Large: 44"-48" hips



### MODEL 898

Basic support provides comfort and a wide range of easy adjustability. 10" elastic back panel with flexible boning. Velfoam sides and front, tapers from 10" to 6", with front contact closure. Tapered elastic side pulls with hook and loop attachments.

Sizes - One Size fits 32"-44" hips



### MODEL 784-HS

Sheer, lightweight form fitting design specially cut for the female figure. White two-ply power net elastic, tapering from 13.5" in back to 7.5" in front - cut with 8" hip development. One pair of 11" removable steels in back. Flexible boning throughout to prevent rolling.

Sizes - Small: 32"-35" hips  
Medium: 36"-39" hips  
Large: 40"-43" hips  
X-Large: 44"-47" hips



### MODEL 896-HS

Combination of biomechanical design and materials provide for comfort and effective support. Heavy duty knitted elastic, tapering from 10" in back to 6.5" in front. One pair of 8" removable stays in back. V-shape heavy duty elastic support pulls with contact closure cinch strap in front. Flexible boning throughout.

Sizes - Small: 30"-32" hips  
Medium: 34"-36" hips  
Large: 38"-40" hips  
X-Large: 42"-44" hips  
2X-Large: 46"-48" hips

# ELASTIC LUMBOSACRAL, DORSOLUMBAR SUPPORTS



## MODEL 2868



Lightweight and comfortable, easy to apply and remove, a simple design and very cost effective. 10" elastic back tapers to 9" in front. 3" criss-cross pattern in back coupled with side pulls provides for three way compression over injured soft tissues. Hook and loop fastening in front and on side pulls; flexible boning throughout.

- Sizes – Small: 30"-34" hips  
 Medium: 35"-38" hips  
 Large: 39"-42" hips  
 X-Large: 43"-46" hips



## MODEL 2869



Same materials, design and application as 2868 above, but cut for a shorter figure. 8" back tapers to 6.5" in front.

- Sizes – Small: 30"-34" hips  
 Medium: 35"-38" hips  
 Large: 39"-42" hips  
 X-Large: 43"-46" hips



## MODEL 889-HS

Cut for the male figure; especially recommended for those who may have trouble applying the more complicated side lace corset type support. White power net elastic, tapering from 20" in back to 9" in front. Contact closure in front. One pair of 18" removable steels in back. Flexible boning throughout garment. Detachable and adjustable criss-cross straps with truss buckle fastening.


- Sizes – Small: 30"-34" hips  
 Medium: 34"-38" hips  
 Large: 38"-42" hips  
 X-Large: 42"-46" hips



## MODEL 777-HS

Same general design and features as 889-HS, but fashioned with 6" hip development for the female figure. Excellent support for those suffering from osteoporosis of the thoracolumbar spine. White power net elastic, tapering from 18" in back to 7" in front. Contact closure in front. One pair of 16" removable steels in back. Flexible boning throughout garment. Detachable and adjustable criss-cross straps with truss buckle fastening.

- Sizes – X-Small: 28"-31" hips  
 Small: 32"-35" hips  
 Medium: 36"-39" hips  
 Large: 40"-44" hips  
 X-Large: 45"-48" hips

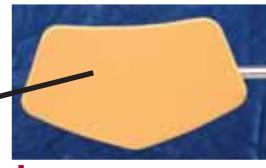

## MODEL 885-X

Provides for total contact support – a gentle well balanced pressure over the lumbar spine.

Simple to fit and easy to apply: a rigid thermoplastic pad is heated to pliability, then molded to the specific contours of the wearer's back. Once the pad is set it is inserted into the pocket of a comfortable elastic support belt. Garment features hook and loop fastening panel.

- Sizes – Small: 24"-31" hips  
 Medium: 28"-35" hips  
 Med. Lge.: 32"-39" hips  
 Large: 38"-45" hips  
 X-Large: 44"-51" hips

